

VEMO

Continuous ECG Recorder


- Quick Guide -

This is veterinary use only.
Please refer to the operation manual before use.


vemo@ebionet.com
Bionet Co., Ltd. - 5F, Dreammark1 Data Center 61 Digital-ro 31 gil, Guro-gu, SEOUL, 08375
+82-2-6300-6410

support@bionetus.com
Bionet America Inc. - 2691 Dow Ave. Suite B, Tustin, CA 92780
+1-877-924-6638


VEMO System component 1




Charging cable



Charging cradle





VEMO ECG sensor




Hydrogel Adhesive 2packs
(1pack / 4pcs)

App store : Bionet VEMO

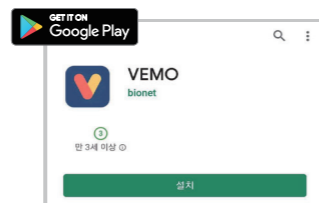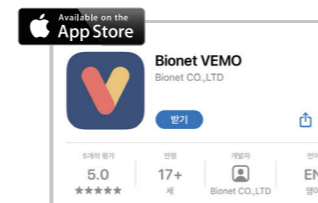

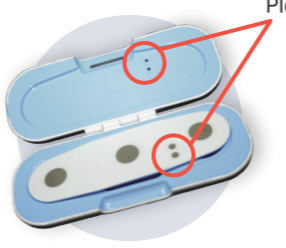

Charging VEMO 2




Please check!


1. Make sure two charging points are connected.



2. Place the VEMO sensor in the charging cradle.



3. Close the cradle.




4. Connect the USB cable to charging cradle.

- Light blue : charging
- Light green : fully charged


How to Authenticate VEMO 3

How to authenticate VEMO


When downloading the VEMO app from the app store, it is set up to process authentication. To use the VEMO app, the wireless, wearable sensor is used for authentication.




1. After downloading the app, please input the hospital's information.



2. Place the wireless sensor on the tablet's screen for the authentication process.



3. After authentication is completed, go to the main screen like the picture on the left.




4. If the authentication was not completed, the app cannot be used & an error message will show up like the picture on the left. (press the back button to restart the authentication process)


Pairing VEMO 4

Steps for pairing


* Caution - Several VEMO sensors around the tablet or mobile can cause pairing error.

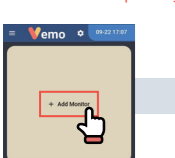


1. Press +Add.

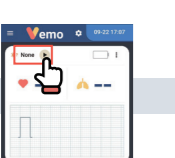


2. Press play button.

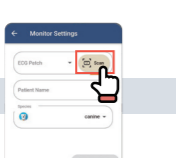


3. Press 'Scan' button.




4. Heart icon window will pop up.

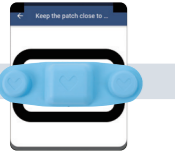

5. Place the VEMO sensor on the heart icon window.

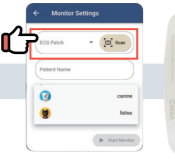


6. Check if the sensor serial number is registered successfully.
Ex) C740139

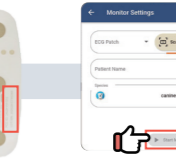


7. Select the type of animal and 'Start Monitor' button.


When 'Auto Scan' has a problem (manual pairing)



1. Check the sensor's serial number.

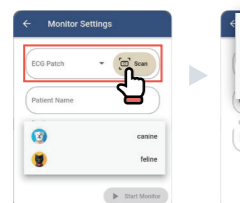


2. Press 'scan' button.




3. Select the serial number.

Applying VEMO 5


Steps for applying adhesive to VEMO sensor

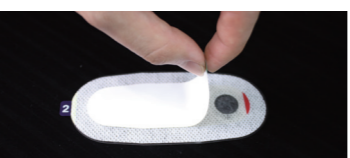


1. Peel out adhesive No.1.



2. Attach the VEMO sensor on the No.1 side.




3. Peel out adhesive No.2.







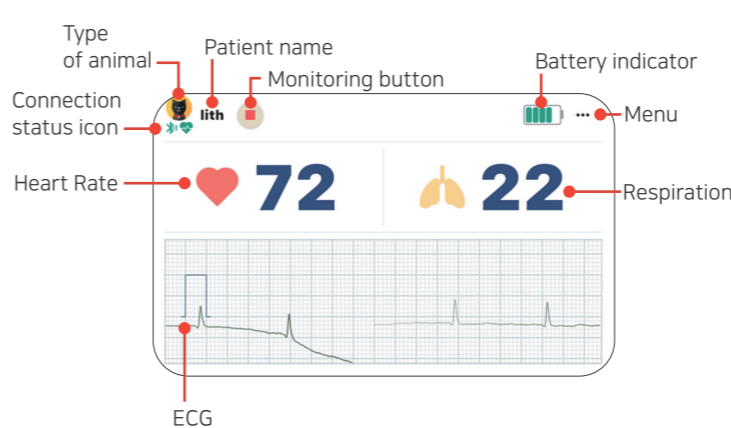

4. Placed the VEMO sensor on the heart with its heart icon 90 degree vertical and clockwise.

Start monitoring 6



Status icon

- Pairing status
 -  Connected
 -  Disconnected
 -  Standby
 -  Error
- ECG lead status
 -  Connected
 -  Error or sensor not applied

Reviewing Recorded Chart & Sharing 7

Chart review

Press Menu and select Vital Signs Chart.







Table chart




Graphic chart

Sharing Chart


Press Share button.

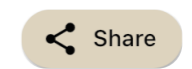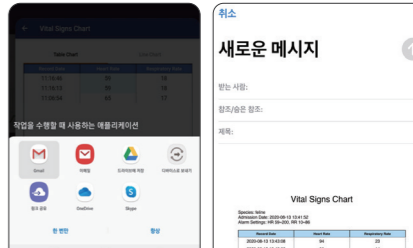
Select sharing method to convert the chart into PDF format.

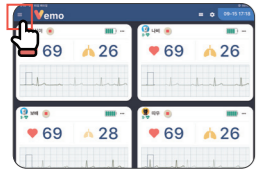



1. Press (☰) button on the top left.

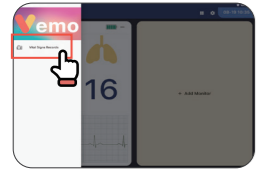


2. Select 'Vital Signs Records' button.

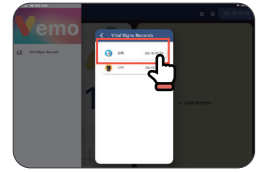


3. Select patient from the list.

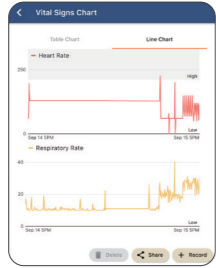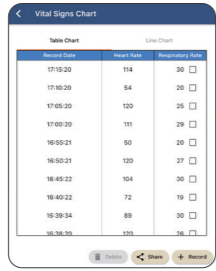

♥ How to set up the alarm

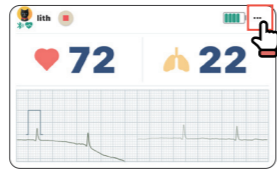


1. Press ... button to select Alarm Settings menu.

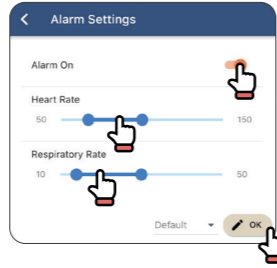


2. Choose Alarm On/Off.
3. Set up the alarm range of Heart Rate, Respiration rate.
4. Alarm will be activated when parameter values are out of alarm range.
5. Press Default setting to return to basic setting.
6. Press OK to save the setting.

♥ Main configuration

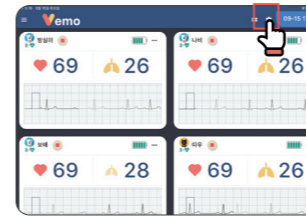


Press the configuration button to see 'Preference'. (⚙)

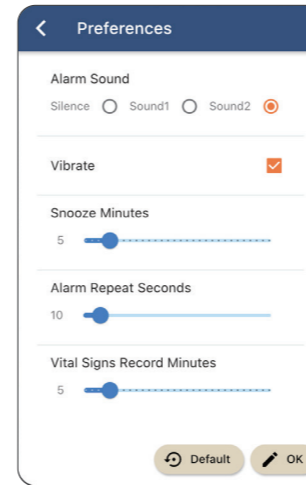


- Alarm Sound: Alarm sound settings, Silence, Sound1/2
- Vibrate: Vibrate On/Off
- Snooze Minutes: Set up Snooze Minutes
- Alarm Repeat Seconds: Set up the alarm notice interval
- Vital Signs Record Minutes: Set up the saving data interval for the Heart Rate, Respiration data


♥ Display layout

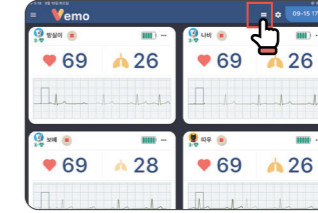

1. Press monitor button at the main display. (📺)

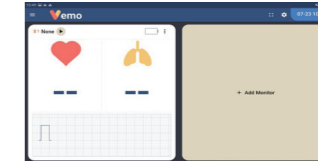

2. By touching the button, display layout will change.

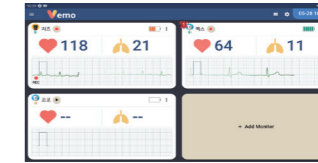

3. Set up the display layout 1X2, 2X2, 2X3 etc.



1. The recommended pairing range between VEMO sensor and tablet is 3m(10ft) up to 10m(32ft).



2. Use VEMO genuine hydrogel adhesive patch to protect VEMO sensor's ECG lead. In case of use other adhesive, warranty will be voided.



3. VEMO sensor is sensitive to heat, magnetic fields, and metal objects. Make sure VEMO sensor is not located near those materials.



4. Check the sensor serial number is read & registered successfully.



5. When VEMO is in use, check pairing for ECG lead status and battery indicator.



Hello,
I'm VEMO!



We deliver highest quality of connected healthcare solutions through constant innovation and sharing the value with our customers.