

Automated Record Keeping for Bionet Veterinary Patient Monitors

B-Link is a software program that enables you to interface Bionet veterinary monitors, including the BM3Vet, BM3Vet Plus and BM5Vet, directly to Windows based computers. The patient data can be saved as a formatted text file or can be pasted directly into the electronic medical record of the veterinary practice management software you already use.

- Paste results directly into ANY Windows based software
- On screen parameter selection
- Set transfer interval in minutes or seconds
- Parameters can also be displayed as a graph
- Data can be saved in formatted text file



B-Link Ver 1.0.0.11


Client / Pet information

Client: <input type="text" value="JOHN DOE"/>	Pet: <input type="text" value="SPARKY"/>	Sex: <input type="text" value="MN"/>
ID: <input type="text" value="3410"/>	Species: <input type="text" value="CANINE"/>	DOB: <input type="text" value="03/04/2002"/>
Phone: <input type="text" value="555-555-5555"/>	Breed: <input type="text" value="SPITZ"/>	Weight: <input type="text" value="345LBS"/>
Misc: <input type="text" value="Over Weight"/>		

Check list BM5Vet.LIF 30 Seconds [Play] [Stop] [Graph] [Refresh]

Date and Time	HR	SPO2-%	SPO2-R	RESP	NIBP-S	NIBP-M	NIBP-D	TEMP1
3/27/2007 5:02:29 PM	120	98	119	30	145	100	85	100.94
3/27/2007 5:02:49 PM	120	98	119	30	145	100	85	100.94
3/27/2007 5:03:19 PM	120	98	119	30	145	100	85	100.94
3/27/2007 5:03:49 PM	120	98	119	30	145	100	85	100.94
3/27/2007 5:04:19 PM	120	98	119	30	145	100	85	100.94
3/27/2007 5:04:49 PM	120	98	119	30	145	100	85	100.94
3/27/2007 5:05:19 PM	120	98	119	30	145	100	85	100.94
3/27/2007 5:05:49 PM	120	98	119	30	145	100	85	100.94
3/27/2007 5:06:19 PM	120	98	119	30	145	100	85	100.94
3/27/2007 5:06:49 PM	120	98	119	30	145	100	85	100.94

Check list chart

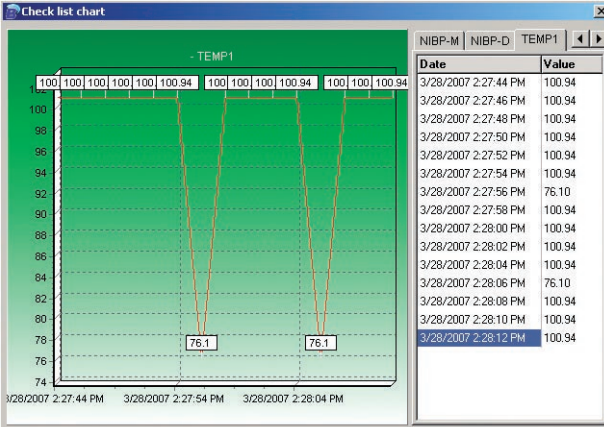

Date	Value
3/28/2007 2:27:44 PM	100.94
3/28/2007 2:27:46 PM	100.94
3/28/2007 2:27:48 PM	100.94
3/28/2007 2:27:50 PM	100.94
3/28/2007 2:27:52 PM	100.94
3/28/2007 2:27:54 PM	100.94
3/28/2007 2:27:56 PM	76.10
3/28/2007 2:27:58 PM	100.94
3/28/2007 2:28:00 PM	100.94
3/28/2007 2:28:02 PM	100.94
3/28/2007 2:28:04 PM	100.94
3/28/2007 2:28:06 PM	76.10
3/28/2007 2:28:08 PM	100.94
3/28/2007 2:28:10 PM	100.94
3/28/2007 2:28:12 PM	100.94

OK Cancel


B-Link Connection Methods



Crossover Cable




Computer to Computer

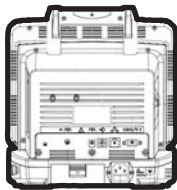

Ethernet LAN Cable



Standard Networking

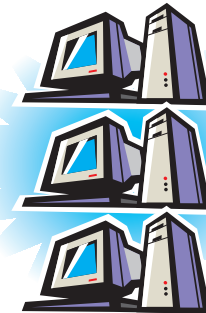

B-Link Parameter Table : Data Transfer by Unit

Parameter	BM3Vet	BM3Vet Plus	BM5Vet	BM5Vet-CS(EtCO2)
HR	●	●	●	●
SPO2-%	●	●	●	●
SPO2-R	●	●	●	●
RESP	●	●	●	●
NIBP-S	●	●	●	●
NIBP-M	●	●	●	●
NIBP-D	●	●	●	●
TEMP1*	●	●	●	●
TEMP2*			●	●
IBP1-S*			●	●
IBP1-M*			●	●
IBP1-D*			●	●
IBP1-P*			●	●
IBP2-S*			●	●
IBP2-M*			●	●
IBP2-D*			●	●
IBP2-P*			●	●
ETCO2				●
ETCO2-R				●
FIO2				●

* Temp and IBP require separate purchase of probe. Please refer to BM3Vet or BM5Vet equipment brochure.